sales@holmarc.com (mailto:sales@holmarc.com) / \* (index.php) \* (mailto:sales@holmarc.com)

(contact\_us.php) / INDIA



Products (products.php) / Physics Lab Instruments (physics lab instruments.php) / Polarization (polarization.php) / Brewster's Angle Apparatus (brewster\_angle.php)






Holmarc's Brewster's Angle Apparatus (Model No: HO-ED-P-01) is designed to study the Brewster's angle phenomenon and the polarization of reflected light. The essential elements of the apparatus consists of a goniometer, a laser light source capable of projecting a light beam that is linearly polarized in its plane of incidence, and a pinhole photo detector with output measurement unit for detecting and measuring the intensity of light reflected. The diode laser and polarizer rotator are mounted on an optical rail.

More

## **Experiment Examples**


- ▶ To measure and plot the graph-reflectivity versus angle of incidence
- ► To find the Brewster's angle (also known as the polarization angle) of glass plate and determination of refractive index 

  +

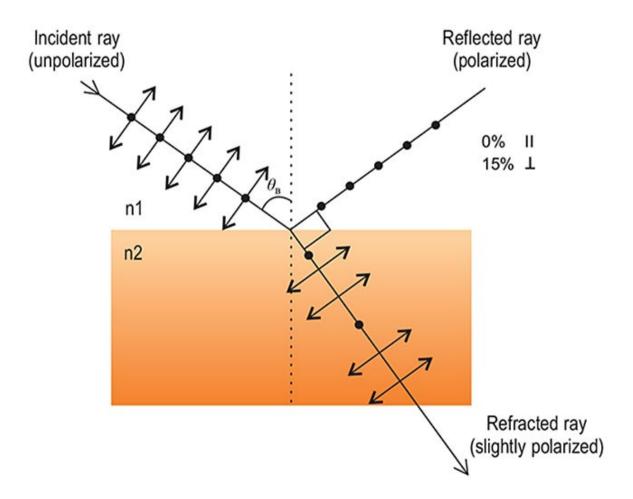

## Related Topics

- Snell's law
- Polarization by reflection
- Fresnel's law of reflection
- Brewster's law

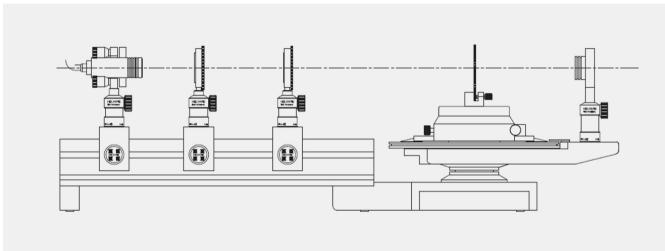

**Mechanical Drawing** 

## Scope of Supply

▶ Goniometer with Detector Mount

+

•	Optical Rail	Т
•	Kinematic Laser Mount	+
•	Polarizer Rotator with Mount	+
•	Glass Slide	+
•	Diode Laser with Power supply (Red)	+
•	Detector Output Measurement Unit	+
•	Accessories	+

(http://www.paypal.co.in/in)

Home (index.php) | Products (products.php) | Supports (supports.php) | Company (company.php) | Contact Us (contact\_us.php)

Privacy Policy | Terms & Conditions

©2016 Holmarc Opto-Mechatronics (P) Ltd. All rights reserved.

Follow us (http://www.facebook.com/Holmarc) (https://twitter.com/Holmarc\_Kochi)

You (http://www.youtube.com/holmarcproducts)